


BEAUTIFIL-Bulk *Restorative*

Advanced Giomer Restorative



CASE STUDY

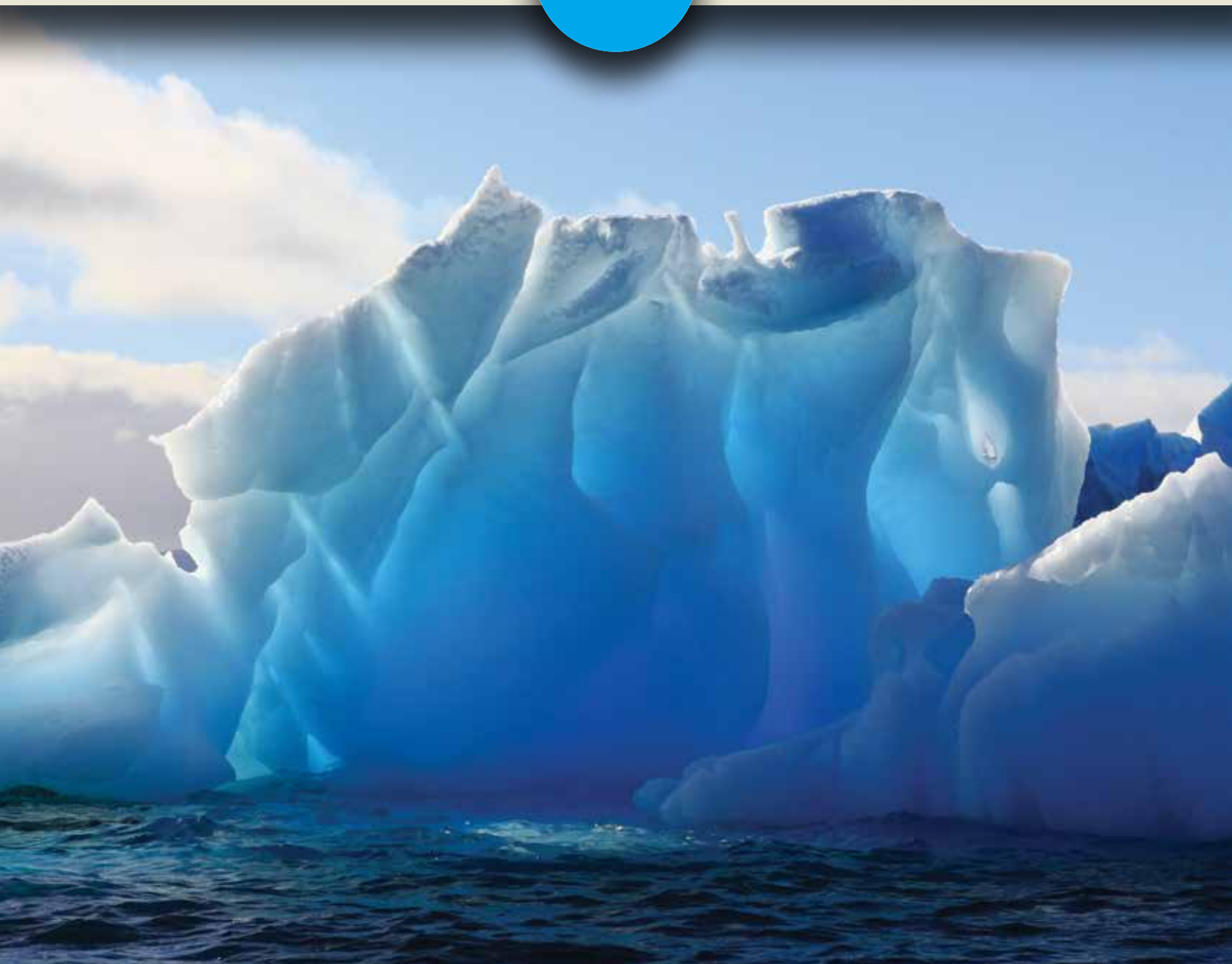

One: Before



Two: After

Photos courtesy of Jack Griffin DDS.

Beautiful-Bulk Flowable incorporates the unique, bioactive benefits of Giomer, a surface pre-reacted glass filler backed by 13 years of clinical success*, within a fast and easy bulk fill composite. Beautiful-Bulk Flowable is designed to overcome the challenges of current bulk fill products, with nice aesthetics, lower shrinkage stress and full polymerization.

*Gordan, V. et al, ed. IADR General Session. Seattle, Washington State Convention Center. March 23 2013. Clinical Evaluation of a Giomer Restorative System: Thirteen-Year Recall. Poster Session.


BEAUTIFIL-Bulk Restorative

Advanced Giomer Restorative

Ideal for posterior restorations, Beautiful-Bulk Restorative is a bioactive packable composite resin that provides exceptional condensability and sculptability. Delivering superb shade stability, low shrinkage stress and full polymerization at 4mm, Beautiful-Bulk Restorative is an aesthetic, advanced Giomer restorative you can rely on.

- High fill ratio (87.0wt% 74.5vol%) reduces polymerization shrinkage (1.7%) and shrinkage stress while increasing compressive and flexural strength.
- Ideal balance of light diffusion and transmission ensures complete polymerization at 4mm depth of cure
- Outstanding Vickers hardness value
- Optimum translucency creates esthetic shades unaffected by surrounding intraoral color
- Sustained fluoride release and rechargability
- High level of radiopacity

Visit www.shofu.com | Call **800.827.4638**

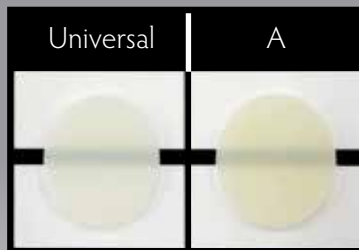

Order Information	
Syringe, 4.5g	
Shades	Item No.
Universal	2034
A	2035
Tips, 0.25g x 20	
Shades	Item No.
Universal	2032
A	2033



Shofu Dental Corporation
1225 Stone Drive
San Marcos, CA
92078-4059, USA

SN202-0718

Advanced Giomer Restorative