BeautiBond Xtreme

Single-Component, Total Universal Bond with Quad-Adhesive Technology



BeautiBond Xtreme

Single-Component, Total Universal Bond with Quad-Adhesive Technology


SERVE PATIENTS INCREASE PRODUCTION GO HOME ON TIME!


Shofu is dedicated to manufacturing the highest quality products utilizing our innovative technologies to enhance minimally invasive cosmetic dentistry.

s h o f u . c o m 8 0 0 . 8 2 7 . 4 6 3 8 customer-service@shofu.com


STRONG DURABLE BOND

Due to **BeautiBond Xtreme's** unique Quad-Adhesive Technology, a single layer gives a strong, fast, and reliable bond. Ideal for direct and indirect restorations using any light-cure, dual-cure or self-cure material without the need for additional primers. **BeautiBond Xtreme** is the truly universal adhesive that offers flexibility with self-etch, selective-etch, or total-etch—the choice is yours!

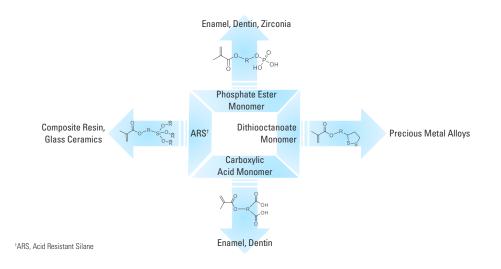


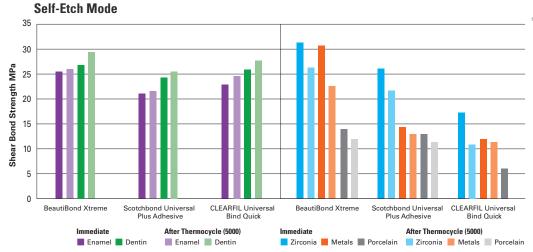

BeautiBond Xtreme is HEMA-free, eliminating the risk of sensitivity and bond hydrolysis. This also improves the stability of the adhesive layer, resulting in a uniform layer that is only 5µm thick and with minimized technique sensitivity. Even after the thermocycling process, no gaps were found between the enamel and dentin adhesive interfaces, allowing for a strong bond and excellent marginal aesthetics. Composite Resin BeautiBond Xtreme Dentin Surface

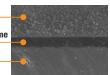

Composite Resin BeautiBond Xtreme Enamel Surface

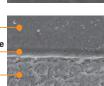



[§]University of Paris, France



PREVENTS SENSITIVITY

BeautiBond Xtreme provides an optimal seal to wet or dry dentin, penetrating deep into the dentinal tubules to create resin tags for secure adhesion, eliminating the risk of postoperative sensitivity and microleakage.

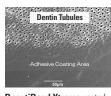

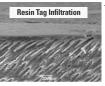

BeautiBond Xtreme coated dentin surface with a single layer application.

Cross sectional SEM image of sealed layer.

AESTHETIC


With extremely low film thickness (5µm), BeautiBond Xtreme eliminates marginal stain lines caused by brands with thicker bonding layers. BeautiBond Xtreme is the best choice when restorations must be picture-perfect.



FAST AND EASY TO USE

BeautiBond Xtreme's ease of use and simplified one-step application makes the bonding process consistently reliable and clinically efficient.

DIRECT RESTORATIONS

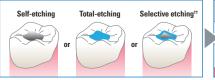

Preparation of the tooth surface Follow the standard method for acid etch. Then rinse and dry

INDIRECT RESTORATIONS



When etching with phosphoric



seconds then blow

with greater force



LED Light: 5 seconds Halogen Light: 10 seconds

until evaporated. acid, rub for 20 seconds.



Apply resin cement, seat and cure as per manufacturer's instruction.

BeautiBond[®] **Xtreme**

shofu.com 800.827.4638 customer-service@shofu.com


RESTORATIVES



BeautiBond Xtreme Unit Dose Set (Item No. Y2702) • (50) 0.1mL unit doses

50 microbrushes




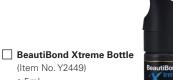



*Shofu internal testing

 Roughen the intaglio restoration surface. Rinse and dry.

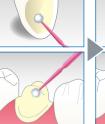



Gently air dry for 3 seconds then blow with greater force until evaporated.

• 5mL

Light-cure. LED Light: 5 seconds Halogen Light: 10 seconds