



# PURE AESTHETICS



UNIVERSAL INVESTMENT



**SHOFU DENTAL CORPORATION**  
Call for Information 800 827 4638 [www.shofu.com](http://www.shofu.com)  
S A N M A R C O S C A L I F O R N I A


## UNIVERSAL INVESTMENT FOR PRESSABLE CERAMICS AND CAST ALLOYS

Ceravety Press & Cast® is extremely versatile and produces perfect fitting casts every time without the need for internal adjustments.

Ideal for traditional castings of metal substructures, pressed porcelain restorations and any alloy from non-precious to high noble. Additionally, the unique stress-releasing formula can be used for pressed to zirconia restorations.

Ceravety Press & Cast mixes easily and has a long working time. A superb universal investment that fulfills all of a technician's requirements and produces clean, precise margins.

Technicians have the freedom to choose between the rapid burn out method and the traditional overnight burn out method. Formulated to reduce stress in pressed ceramics, it can be easily divested.

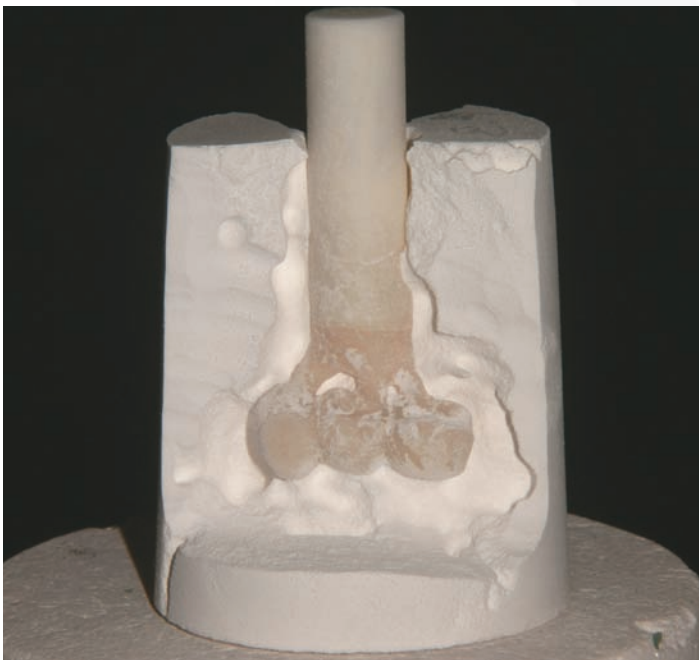

## PHYSICAL PROPERTIES

Powder: Liquid Ratio	100g : 20ml
Working Time (min.)	6.0
Setting Time (min.)	9.0
Thermal Expansion 850°C (1562°F) %	0.9
Compressive Strength (MPa) after two hours	10.0

## CERAVETY – DESIGNED TO BE WORRY-FREE!

- Ceravety Press & Cast** works with all alloys and pressable ceramics
- Ceravety Press & Cast** provides a long working time; use ring or ringless techniques with confidence
- Ceravety Press & Cast** has a quick or overnight burnout technique for full control of casting
- Ceravety Press & Cast** ensures effortless removal and reduces possibility of chipping of ceramic
- Ceravety Press & Cast** has less burning (reaction layer) of pressable ceramics and investment materials
- Ceravety Press & Cast** increases precision of the fit, resulting in a time saving predictable result



### POWDER (PN 6966)

120 (100g) packages

### LIQUID (PN 6967)

2 liters of the Press & Cast liquid



### MANUFACTURER

**SHOFU INC.** 11 Kamitakamatsu-cho, Fukuine, Higashiyama-ku, Kyoto 605-0983, Japan

**SHOFU INC. Singapore Branch** 10 Science Park Road #03-12, The Alpha, Science Park II, Singapore 117684

**SHOFU DENTAL CORPORATION** 1225 Stone Drive, San Marcos, CA 92078-4059, USA

**SHOFU DENTAL PRODUCTS LTD.** Duke's Factory, Chiddingstone Causeway, Tonbridge, Kent, TN11 8 JU, UK

**SHOFU DENTAL SUPPLIES (Shanghai) CO. LTD.** No.645 Jiye Road, Sheshan Industry Park, Songjiang 201602, Shanghai, China

**EU REPRESENTATIVE: SHOFU DENTAL GmbH** Am Brüll 17, 40878 Ratingen, Germany

**SHOFU**



**CERAMICS**